


STANDING OVATION

IS

ZEST FOR LIFE

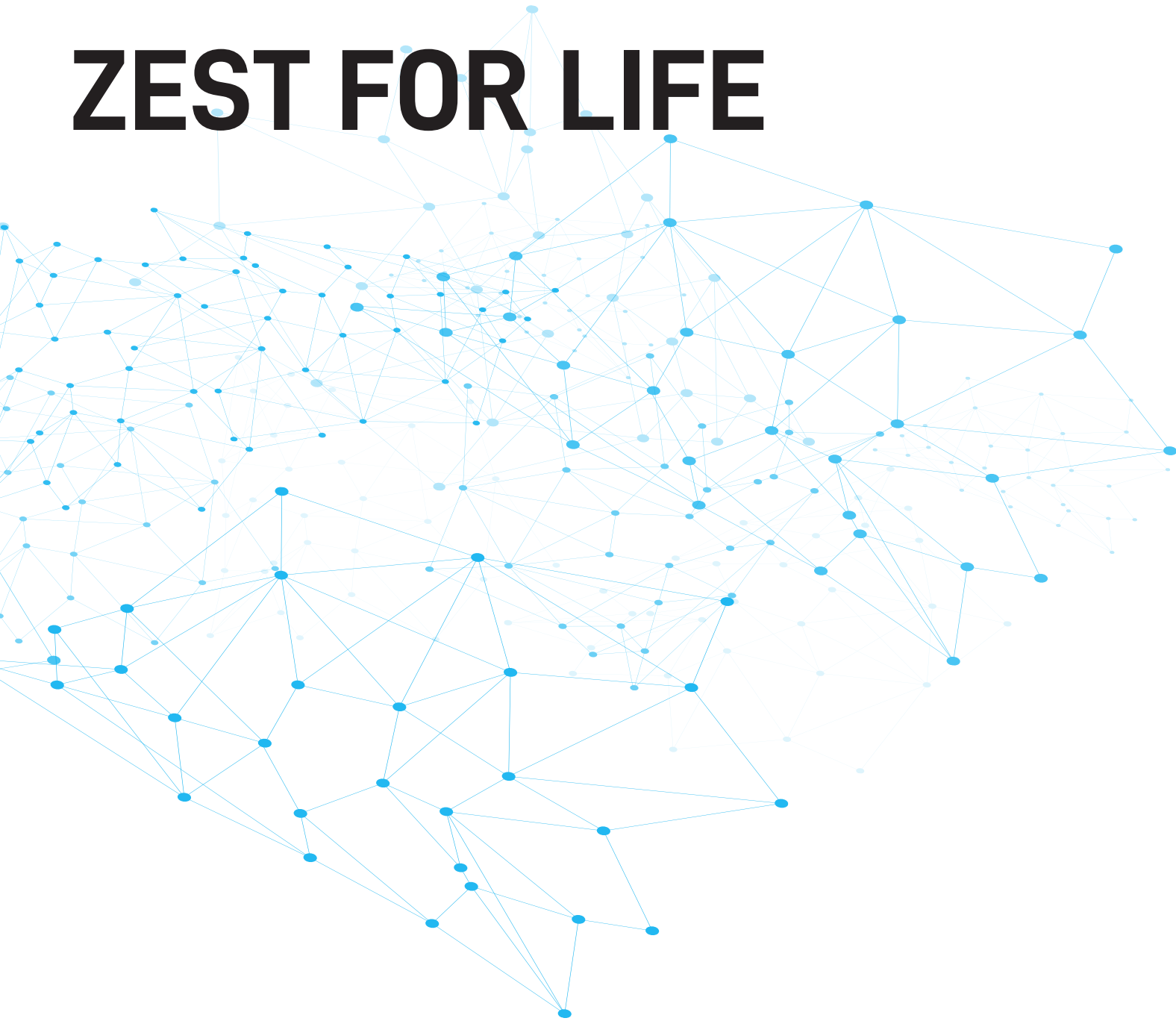

STANDING
OVATION

DEVELOPED FOR & WITH PEOPLE AFFECTED


About ten years ago **Peter Lammer** had a serious motorcycle accident on his way to work, which brought the four-time family father and restaurateur to the brink of despair.

The injuries he suffered affected the strength of his legs severely, resulting in a 80% disability rating. Under these circumstances working in a standing position on a regular basis was unthinkable.



The love for his family as well as the passion for his job were the necessary motivation for the 52-year-old to develop a unique solution together with his partner and friend Bernhard Tichy.

Contrary to all medical predictions, Peter Lammer is now back at work in his kitchen every day. Thanks to **STANDING OVATION**, he can provide for his family, employ his staff and his **ZEST FOR LIFE** has returned.

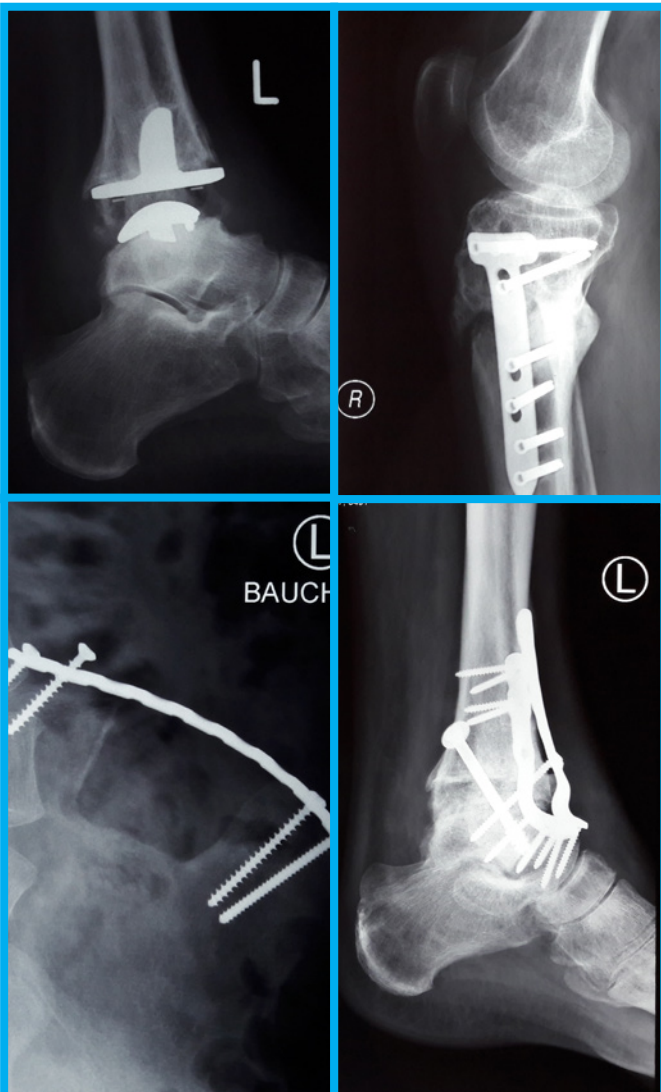

WHAT IS STANDING OVATION?

STANDING OVATION is a **stand and mobility aid** that helps relieve strain on the lower extremities. It consists of three components:

1. rail system

2. C-bar

3. seat lifting unit

The 3-way rail system is either overhead-mounted or a free-standing structure with a C-shaped seat bar and a seat lifting unit attached to it, thereby allowing a three-dimensional movement pattern.

The user's torso weight is transferred to **STANDING OVATION** via the device and thus your hands remain free for work. Likewise, when lifting loads, any additional weight is taken over by the device and there is no additional strain on the legs.

The spring-loaded rotation unit provides the user with maximum freedom of movement and maximum accessibility at the workplace or wherever **STANDING OVATION** is used.

The advantages as demonstrated by Peter Lammer

- lower extremity relief
- less pain, less leg swelling
- reduction of medication
- allows independent picking up and moving of objects (independent work)
- enhances mobility
- increases your work speed
- economic independence
- financial security through independence
- meaningful employment
- **ZEST FOR LIFE**




UNIQUE & SIMPLE TECHNOLOGY..

The rail system


It consists of two longitudinal rails "L" and one rail bridge "B"; curves and circle tracks can be integrated into the construction. The rails are equipped with rolling apparatuses driven by muscle power.

- individual planning & conception
- rail length L is freely selectable
- rail bridge width B up to max. 6 m
- smooth-running rolling apparatuses

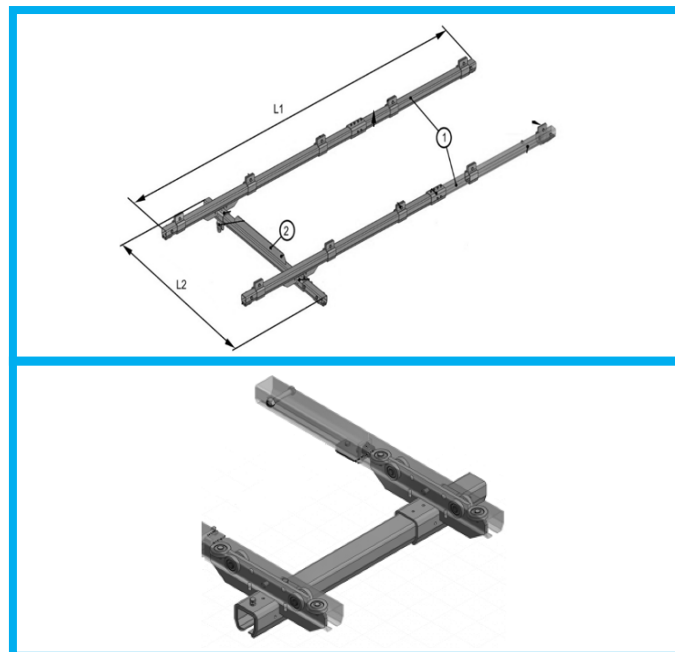

The C-bar

- ergonomically shaped and removable
- individually configurable
- unrestricted axial rotation
- high seating comfort with loose-control
- fall arrester
- standard load capacity up to 125 kg (special designs also possible)

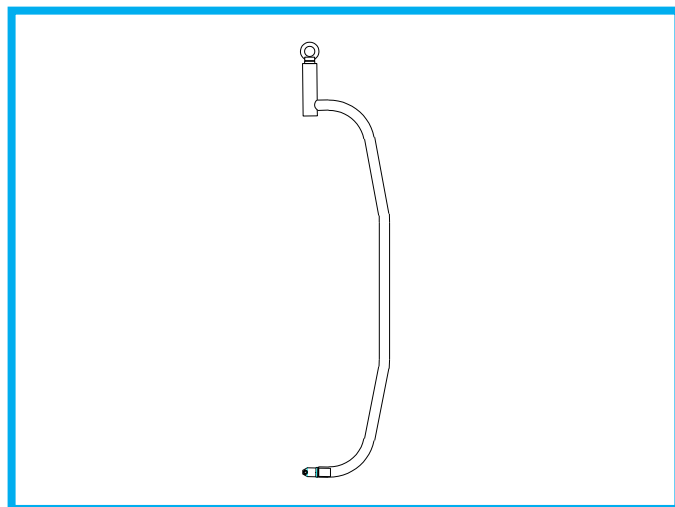

The seat lifting unit

- seat post with gas pressure spring
- one-hand operation
- freely adjustable relief of the legs
- stepless lifting and lowering
- various saddle solutions, also special medical saddles




..THIS IS **STANDING OVATION**

The freestanding carrying frame:

Thanks to the carrying frame, any place of work can be adapted even without a suitable ceiling height or building fabric.


- individual planning & conception
- rail length L is freely selectable
- rail bridge width B up to max. 6 m
- height H – 2.4 metres
- simple assembly
- change of the building substance not necessary
- outdoor use is possible

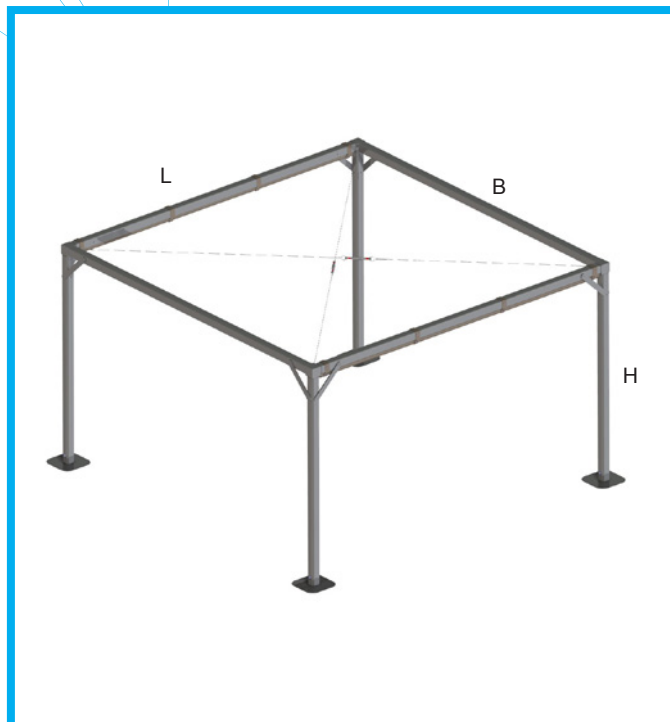

MULTI-ROOM SOLUTIONS ARE ALSO FEASIBLE

The advantages in the workplace


STANDING OVATION will generate added value for employees and employers in the industrial environment and at the workplace. By simply relieving the strain on the legs, many employees can continue or resume their work after injury or illness. It also helps counteract premature signs of wear and tear.

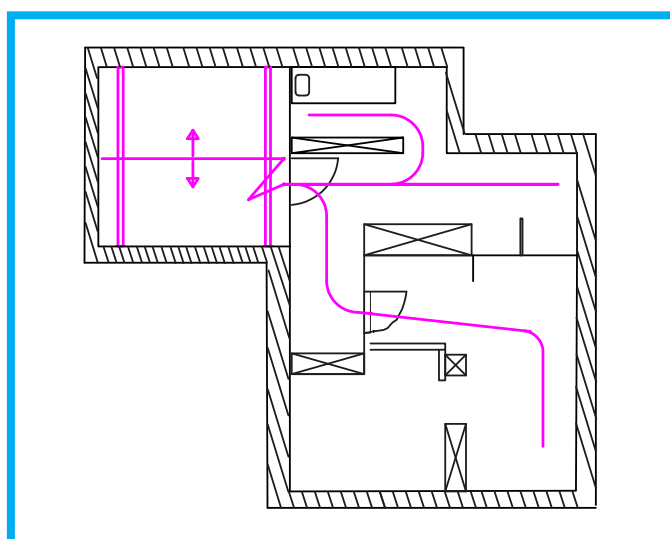

Advantages in rehabilitation

In a therapeutic environment, **STANDING OVATION** (incl. safety fall arrester) will be used as a therapy device to minimise the risk of falling. Postoperatively, **STANDING OVATION** can also facilitate initial mobilisation, e.g. instead of water therapy. Patients and therapists at the AUVA Reha-Institut Bad Häring already have a device at their disposal.

TALK TO US ABOUT ZEST FOR LIFE

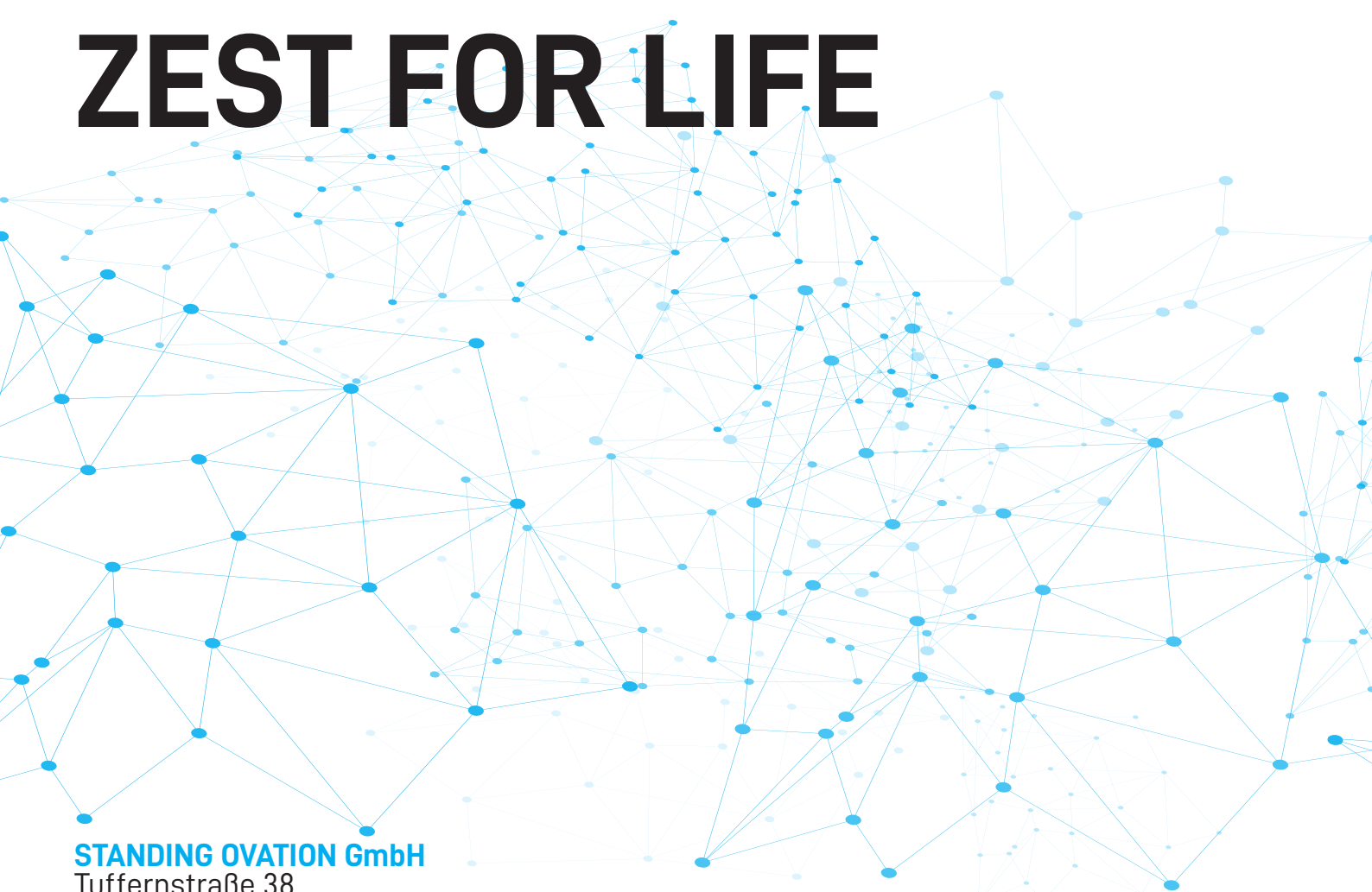

STANDING OVATION GmbH

Tuffernstraße 38
5300 Hallwang

ATU74194123
FN: 507230b

office@standingovation.at
www.standingovation.at

Bernhard Tichy: Information, Technology and Sales
+43 676 70 47 555



STANDING OVATION is patented and protected by trademark law.

Take a look at STANDING OVATION!
Simply scan the QR code and see how Standing Ovation works in practice.
Pat.: # AT519729